


DONJOY®

Performance & Lifestyle

X-Act ROM Knee

Precise Protection.

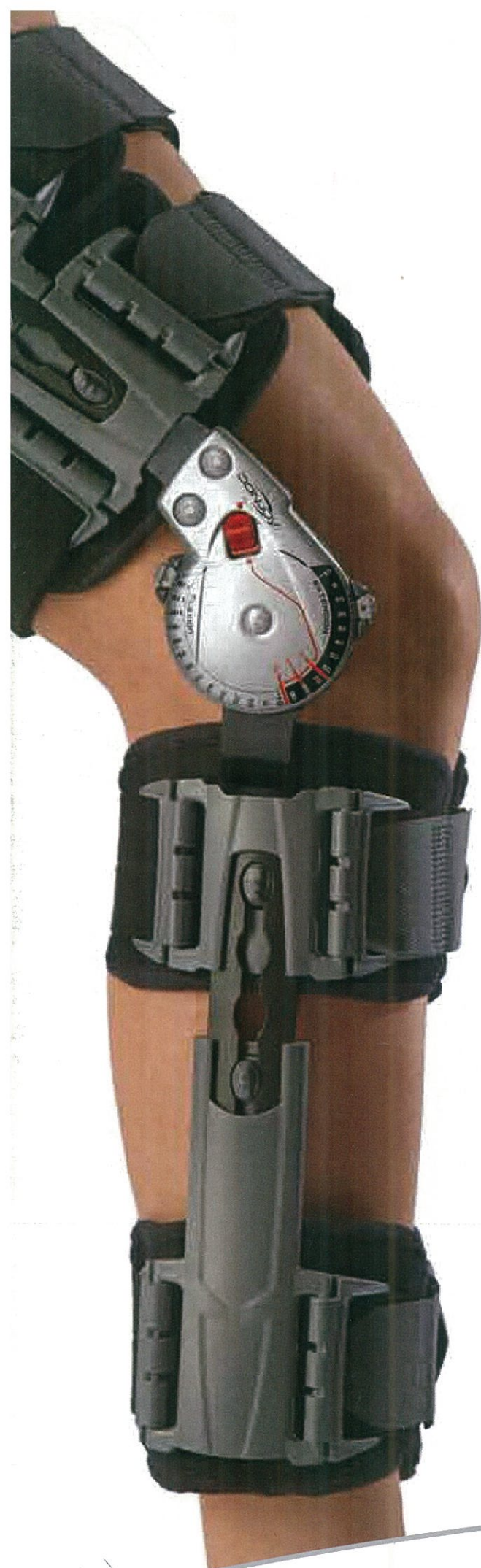

DonJoy's new X-Act ROM Knee Post-Op brace takes comfort, fit and functionality to the next level, providing a better healing environment after surgery.

Indications


Immobilization and protected range of motion associated with ACL, PCL, LCL, and MCL surgeries and meniscal repairs

Durable Aluminium Hinges for Precise Protection

- The new, streamlined, durable aluminium hinges provide precise protection and range of motion control in 10° increments
 - Extension settings from -10° to 90°
 - Flexion settings from -10° to 120°
 - Easy Quick-lock at 0°, 15° and 30°

New Features

Buckles attach quickly and swivel around the leg for immediate contouring with minimal application time.

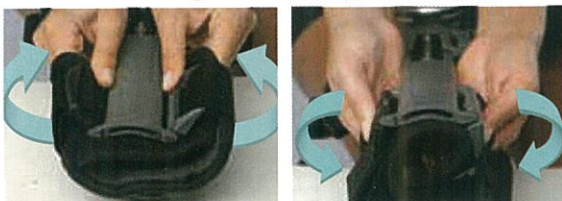

Push clip buckles snap into place for faster fit.



Quick-clip Buckles for Ease of Application

- The swivel action quick-clip buckles attach easily and provide a close-fitting contour around the leg with minimal application time



The X-Act ROM is DonJoy's lightest post-op brace. It is quick and easy to apply, offers a comfortable, contoured fit and provides precise protection during healing.

Tele-Fit Mechanism for Optimal Fit

- Allows independent length adjustment of the proximal and distal straps, making the X-Act Rom Knee a universal fit
- Telescoping adjustability from 48 cm to 63.5 cm (19" to 25") to accommodate a wide variety of patients' sizes
- All 4 sliding cuffs telescope independently allowing placement of straps away from surgical site

Self-adhesive Foam Liners for Brace Stability

- Unique, self-adhesive foam lining provides sticky contact points between the leg and the brace to reduce slippage





X-Act ROM Knee

ORDERING INFORMATION

PART NUMBER	DESCRIPTION	SIZE
11-2151-9	X-Act ROM Knee Brace	Universal *
ACCESSORIES		
25-6755	Brace Lock Kit	-

Length Adjustability: cm to cm [" to "]

*Maximum Thigh Circumference: cm ["]

DJO Global | a Guildford Business Park | Guildford | Surrey | GU2 8XG | ENGLAND

www.DJOglobal.eu



Exclusively Distributed by:

Syarikat Wellchem Sdn. Bhd.

928-929, Jalan 17/38, 46400 Petaling Jaya, Selangor D.E., Malaysia

Tel (03) 7957 2449/ 7957 0218 info@wellchem.com

www.wellchem.com

Together in Motion™